

***Colgate***<sup>®</sup>

powered by **OMRON**

ProClinical<sup>™</sup> **A1500**

**INSTRUCTION MANUAL**

ProClinical™ **A1500**

INSTRUCTION MANUAL


# ENJOY YOUR BRUSH -

It Automatically Adjusts Speeds and Cleaning Actions for Healthier Teeth and Gums.

Please follow all instructions. Save for future use.

## PARTS

---

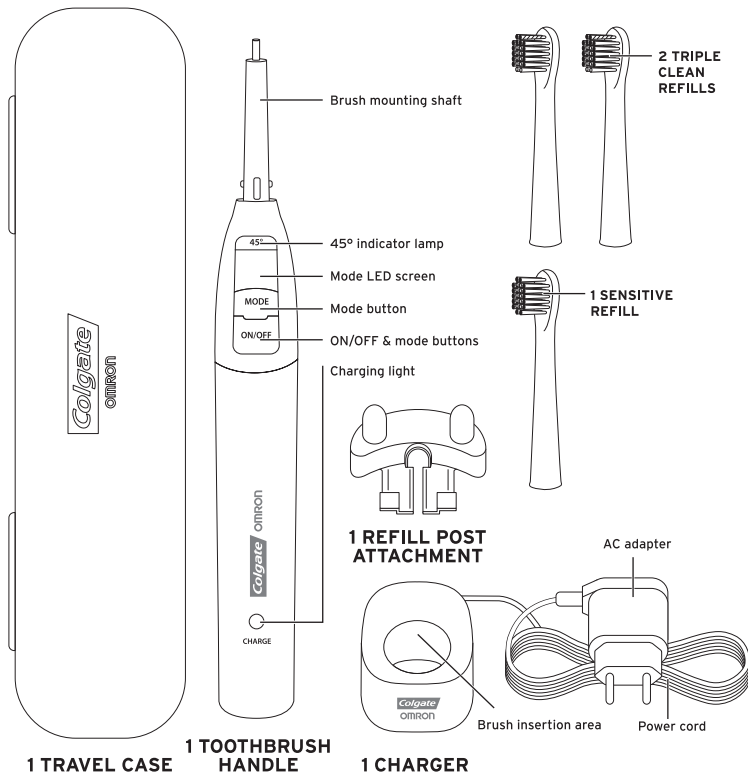

# IMPORTANT SAFEGUARDS

---

## READ ALL INSTRUCTIONS BEFORE USING

Basic safety precautions should always be followed including the following:

### DANGER

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not touch the AC adaptor with wet hands.
- AC adaptors and chargers should always be completely dry when connecting to power.

### WARNING

- The Supply Cord cannot be replaced. If the cord is damaged the appliance should be scrapped.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Use the toothbrush with brush heads and charger provided by the manufacturer only.
- UK (only) supply connection: Type C plug provided is intended for use in a Shaver Socket. It is not suitable for connection to a UK mains socket unless used with an approved (e.g. BSI or ASTA) FUSED Converter Plug with a 3A max fuse or Shaver Adaptor with a 1A fuse.
- Use the charger only within rated capacity of outlets or voltages specified.
- Keep product away from heated surfaces.
- Use this product only for its intended use described in this manual.
- No user serviceable parts. Do not disassemble the handle except for battery disposal, as it may cause personal injury or device failure.
- Stop using this product if it is damaged in any way (cord, plug, brush head, toothbrush handle, charger, and/or AC adaptor).
- Not for children under 3. Keep out of reach of small children.
- If after using this toothbrush, you experience excessive bleeding or if bleeding continues to occur after a week of use, please contact your dentist as it may be an indication of a more serious oral care issue.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.

### NOTICE

- For energy savings purposes, charge only when the battery level is too low. When not in use, keep the charger unplugged from the socket.


### Electromagnetic Compatibility

The Colgate® ProClinical® electric toothbrush complies with all standards regarding electromagnetic fields. It is safe to use if handled properly according to the instructions provided in this manual.

## SAVE THESE INSTRUCTIONS

## 1

## QUICK START

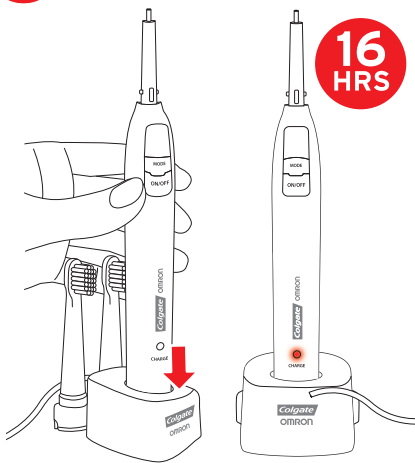


## CHARGE

Brush comes partially charged for initial use. For best results, charge 16 hours before regular use. Adaptor may be required (not included).


- Put handle into charger.
- Plug into outlet for **16 hours** to ensure a full charge.
- **Light** on handle is **red** when charging.

## 2

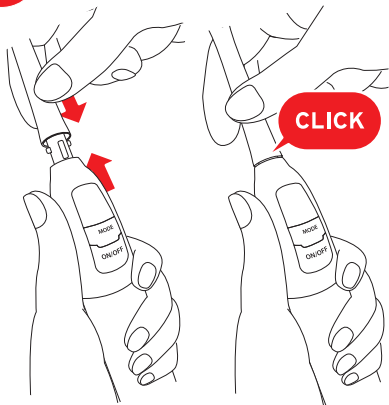




## ATTACH BRUSH HEAD


- **Slide brush head** onto handle.
- Match **grooves**.
- Push until you feel a **click**.

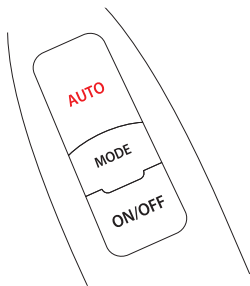

## **AUTO MODE - 3 SPEEDS & 3 BRUSH ACTIONS**


- Screen shows **MODES**.
- Press **ON** (AUTO on screen).
- **Rotate** brush to **listen** for **speed** changes.

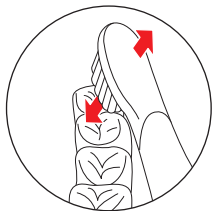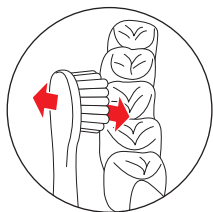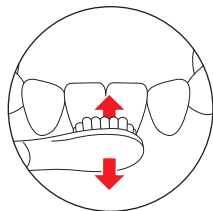

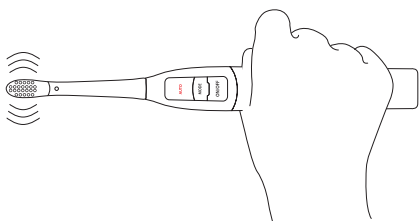


- **Low-speed to gently brush tooth surfaces.**

Brush moves **side-to-side**.

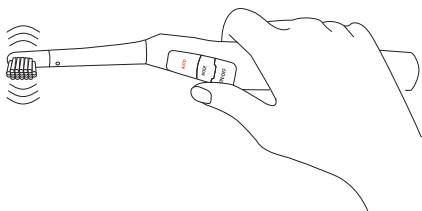

- **Mid-speed when holding handle at 45° to brush at the gum line.**


Brush **pivots** from side-to-side and up & down to clean **teeth and gums**.

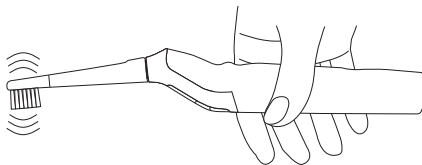

- **High-speed to clean behind & between teeth and molars.**

Brush moves **up & down**.

4

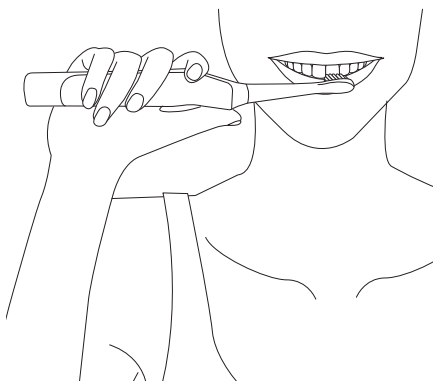

## BRUSH

- **Wet brush, apply toothpaste.** Put brush **in mouth**, press **ON**.
- **Move brush slowly** over teeth. Brush **sensors** work best when **your head is upright**.
- **Timer** alerts every **30 seconds** to brush different area of mouth.
- **Brush stops** after **2 minutes**. **Rinse brush**. Wipe or shake off water.

5



## CHARGE

- **Charge 16 hours** once a week or before if brush is operating slower.
- Leave charger unplugged between charges.



# MODE & SPEED FOR BEST CLEANING EXPERIENCE

---

- **See and Hear 3 MODES: Auto, Optimum, and Deep Clean.**
- **MODES** combine the right **speed** with the best **cleaning action**.
- Press **ON** (MODE last used shows on screen).
- Press **MODE** on handle to select **Auto, Optimum, or Deep Clean**.
- Use thumb indents to help position at 45° (recommended angle to clean along the gum line).
  - Blue light** is ON when 45° position is achieved.

## AUTO MODE (RECOMMENDED)

For new users, easing into your brushing experience in Auto Mode can take up to 2 weeks.

Automatically adjusts speeds and cleaning actions for Healthier teeth and gums.

- **Slower**, side-to-side strokes gently clean **tooth surfaces**.
- **Mid-speed**, up-down and side-to-side strokes clean teeth and gums where plaque builds up most: along the **gum line**.
- **Faster**, up-down strokes clean deep **behind** and **between teeth and molars**.

## OPTIMUM MODE


- **Mid-speed**, up-down and side-to-side strokes clean teeth and gums where plaque builds up most: along the **gum line**.

## DEEP CLEAN MODE

- **Faster**, up-down strokes clean deep **behind** and **between teeth and molars**.

## REFILL BRUSH HEADS


---

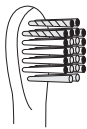

### TRIPLE CLEAN

Cleans in 3 ways

- Spiral bristles help to remove surface stains.
- Tapered outer bristles clean deep along the gum line.
- Inner bristles remove plaque from tooth surfaces.

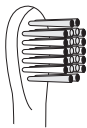

### SENSITIVE

Gentle and Deep Cleaning


- Ultra soft tapered bristles clean gently along the gum line.
- Extra fine inner bristles gently clean tooth surfaces.

Replace with a new Colgate® ProClinical™ Refill brush head **every 3 months**.

## CLEANING BRUSH AND CHARGER

---

- Wash brush head with water.
- Wipe handle and charger with damp cloth.
- Dry brush head and charger completely before putting handle back in charger.
- Never put charger in water.

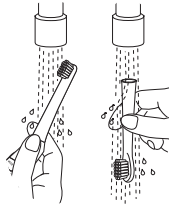

## CHARGER - REFILL ATTACHMENT

---



## BATTERY REMOVAL FOR RECYCLING


---



### REMOVE BATTERY BEFORE DISPOSAL

Remove recyclable rechargeable battery in compliance with environmental regulations before disposing of toothbrush. Do not open handle unless permanently disposing (will void warranty).

### TO REMOVE BATTERY FROM HANDLE:

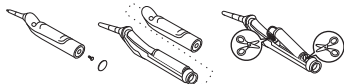

1. Remove cap.
2. Remove screw on bottom of handle with screwdriver.
3. Lift off handle cover.
4. Peel off tape to remove battery. Cut cords.
5. Tape end of battery to prevent shorting.

## CORRECT DISPOSAL



Do not throw away this product with normal household waste at the end of its life. Always remove battery before discarding. Please use a return collection system to ensure responsible and environmentally safe recycling of the device and the battery.

## TROUBLESHOOTING

| PROBLEM   | LIKELY CAUSE   | SOLUTION   |
|---|--|--|
| Tickling sensation while brushing   | It is normal to experience a slight tickling sensation when you first start using this brush.        | After a week you will become more accustomed to this brush and the sensations will diminish. |
| In Auto Mode, the brush seems to change modes automatically for no reason | The brush's Sensor technology works best when head is upright.                                       | Hold your head up right when brushing.   |
| Brush head rattles or is unstable in use                                  | Brush head is not correctly attached to the handle.  | Re-attach brush head by matching grooves and push until you feel a click.                    |
| Slow or no movement   | Battery level is too low.  | Place in charger and charge for 16 hours before use.   |
| Toothbrush and charger are warm   | It is normal for the brush and charger to slightly warm while charging.                              | N/A  |
| Charge indicator light does not light up when brush is placed in charger  | Charger is not properly plugged in or outlet has no power.   | Check that the charger is properly plugged in and socket is live.                            |
| Toothbrush is charged but loses power after a short time of use           | For first time use or use after a long period of not using, the battery requires 16 hours of charge. | Place in charger and charge continuously for 16 hours.                                       |

If your problem is not solved by the suggestions above, please contact your Colgate Consumer Affairs Centre, 00-800-321-321-32.

## 30 DAY CHALLENGE

Try Colgate® ProClinical™ toothbrush for 30 days from date of purchase. If you are not satisfied with this product, return the product within 45 days of purchase for a refund of the purchase price (shipping and handling excluded). You must return the product in its original box with all contents: handle, refill brush heads, charger, travel case (if applicable), together with the original cash register receipt showing this product was purchased in the United Kingdom. Please allow 6-8 weeks for handling. To find the return address closest to you, visit [www.colgateproclinical.co.uk](http://www.colgateproclinical.co.uk). For further details, please contact Colgate Consumer Affairs Centre, 00-800-321-321-32.

Consumers have statutory rights under applicable national laws relating to the sale of consumer products which are in no way affected by the above.

## WARRANTY

Colgate-Palmolive UK Limited grants a 2 year warranty on this product from the date of purchase within the United Kingdom. Within the warranty period, we will repair or replace the product (the choice of whether to repair or replace being ours) in the unlikely event of failure or defects resulting from faults in materials or manufacturing. Following any replacement, the warranty will continue from the date of original purchase only. This warranty does not cover brush heads, or any damage due to improper use (any use not in accordance with the User Manual Instructions and Warnings) or normal wear or use.

To obtain service within the warranty period, please contact Colgate Consumer Affairs Centre, 00-800-321-321-32.



**COLGATE-PALMOLIVE**

© GUILDFORD GU2 8JZ

® DUBLIN 24, IRELAND

[www.colgateproclinical.co.uk](http://www.colgateproclinical.co.uk)

*Made in China under license of  
Colgate-Palmolive Company, Trademark Owners*

*Version 01, January 2012*